

GASTROENTEROLOGY



An Indian company
that is raising the
standards of
medical device
manufacturing
in India

BioradMedisys is a leading global speciality manufacturer of a large portfolio of medical devices (Indovasive) and orthopedic implants (Orthovasive).

Indovasive, a division of BioradMedisys Pvt Ltd, is one of the most promising medical devices manufacturing set-ups in India, offering a complete range of disposables in endourology (URS, PCN, RIRS, TURP procedures), gastroenterology, interventional radiology, gynaecology and vascular access segments.

We are market leaders in surgical disposables in India, and striving to become the market leaders for the world as well.

The company has state-of-art manufacturing facilities (as per GMP and WHO Standards) at Bangalore and Pune, equipped with cutting-edge technologies and sophisticated quality control departments and labs.

We are an ISO 9001:2008, ISO 13485 and CE certified company, and are determined to reach greater heights with new certifications for quality. We have a presence in more than 50+ countries around the globe with our high quality products.



Table contents

Tissue Acquisition

Biopsy Forceps	3
----------------------	---

Hemostasis

Ligasmart Bands	5
Klosure Clips	5
Sclerotherapy Needle	6
Bioqlot Powder & Pump	6

Dilatation

Xpand Rigid Dilator	8
Spring Tip Guidewire	8
Xpand Multistage Balloon Dilator	9

Biliary Devices

Guidewires	11
Sphincterotome	11
Cytology Brush	12
Stone Extraction Balloon	12
Stone Extraction Basket	13
Lithotripsy Basket	14
Lithotripter	14
Tefloflex	14

Self Expanding Metal Stents (SEMS)

Esophageal Partially Covered Stent	16
Esophageal Fully Covered Stent	16
Biliary Uncovered Stent	17
Biliary Fully Covered Stent	17
Biliary Partially Covered Stent	17

Endoscopic Ultra Sound (EUS)

FNA Needle	19
FNB Needle	19

Enteral Feeding

PEG Kit	21
Jejunal Tube	22

3rd Space Endoscopic Accessories

Polypectomy Snare	24
-------------------------	----


Miscellaneous

Loop Basket	26
Bio-NBDC Catheter	26
Bio-Bilbao Dotter	27

TISSUE ACQUISITION




- Dual Pull wire helps in achieving a tangential bite
- Fenestrated cup design facilitates minimum blood interference with tissue sample
- Shaft covered with coated sheet that facilitates smooth passability of forceps through the working channel
- Available in a wide variety of jaw designs

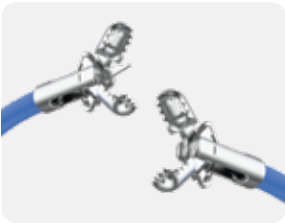


Alligator Cup Forcep

Available In: With Needle, Without Needle
Jaw Diameter: 2.3mm
Working Length: 160cm,230cm,150cm
Jaw Diameter: 1.8mm
Working Length: 150cm

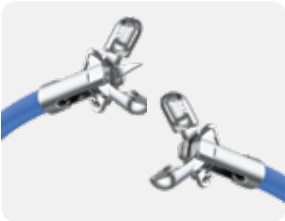


Oval Cup Forcep

Available In: With Needle, Without Needle
Jaw Diameter: 2.3mm
Working length: 160cm, 230cm, 150cm
Jaw Diameter: 1.8mm
Working Length: 150cm

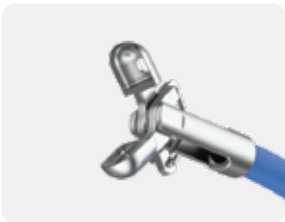

Standard Cup Forcep

Available In: Without Needle
Jaw Diameter: 2.3 mm
Working Length: 160cm, 230cm
Jaw Diameter: 1.8mm
Working Length: 150cm

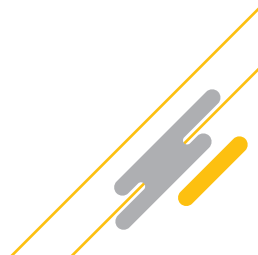

HEMOSTASIS



LIGASMART

Multiband Ligator

- Handle designed to provide an audible and a tactile feedback during the deployment of bands
- Bands made of latex free material to reduces risk of latex allergy
- Each band is square moulded to provide dependable performance



Multiband Barrel

Band Barrel Available In 4, 5, 6 Bands

Multiband Ligator Kit

Band Barrel Available In 4, 5, 6 Bands

Multiband Ligator Kit

Band

Barrel

Handle

Band Barrel Unit

Trigger Thread

Hook Tube.

KLOSURE

Effective And Quick Hemostasis

- Pre-loaded clips provide faster response time in emergency bleeding procedures
- 15mm open jaw diameter helps in good tissue purchase & better hemostasis
- Repositionable design for accurate placement
- 360° rotation helps in deployment of clip in tortuous locations
- MRI compatible clip



Available In:

Size: 15 mm


Lengths: 230cm



TRUJET

Sclerotherapy Needle

- Unique catheter design that helps needle actuation in tortuous anatomy
- Locking mechanism to secure needle advancement during injection
- Transparent outer catheter sheath allows visualization during needle advancement

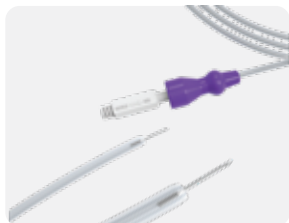


Available In:

Sizes: 21G, 23G, 25G

Lengths: 200cm, 240cm

BIOQLOT

- It's a non thermal method of Hemostasis for non variceal bleeds.
- No tissue trauma as compared to other methods.
- Powder sprayed through a pump, covers a large surface area and hence does not require specific targeting like other modalities.
- Made from sea shells that provide quick hemostasis.

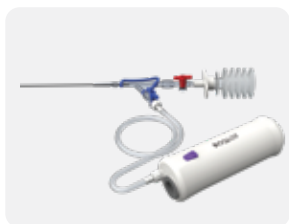

DILATATION



XPAND

Over The Guide Wire
Esophageal Dilator

- Flexible shaft designed for improved precision and control
- Tapered tip with 2 radio opaque markers placed at the distal and proximal end of the tip to enhance its visibility



Type: Bougie Dilator Set

Standard Accessories In Set : Bougie Dilators,
Spring Tip Head Guidewire And Cleaning Brush
Size Available Per Set: 5mm-15mm
Dilator Units Per Set: 7 Units.
Working Length: 70cm, 100cm

Bougie Dilator

Size Available: 5mm, 7mm, 9mm, 11mm,
12.8mm, 14mm, 15mm
Working Length: 70cm, 100cm



Accessories For Dilator

Available: Spring Tip Guidewire
Diameter: 0.038"
Working Length: 250cm

Available: Cleaning Brush


Working Length: 120cm



XPAND

Balloon Dilator

- Single balloon gives 3-distinct dilation diameters at 3 distinct pressures
- Balloon made of non-compliant material for uniform dilation
- 2-Radiopaque markers for X-ray visualization
- Rounded shoulders designed to facilitate visualization of stricture while dilating

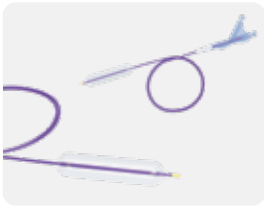

Multistage Over The Wire Balloon

Available In: 7 Fr

Guidewire: 0.035

Shaft Size: 6.8fr

Balloon Diameter: 6-7-8mm, 8-9-10mm, 10-11-12mm,

12-13.5-15mm, 15-16.5-18mm

Balloon Length: 5.5cm

Balloon Catheter Length: 180cm, 240cm

BILIARY DEVICES



CHASE

Guidewire

- It's an access and exchange wire
- Atraumatic tip design aids in easy passability through challenging anatomy and tight strictures
- Engineered with a core alloy providing an excellent fluorovisibility



Available In :

Tip Type: Straight, Angled

Diameter: 0.018", 0.025", 0.032", 0.035"

Wire Length: 260cm, 450cm

TRUSLIT

Sphincterotome

- It's a triple lumen sphincterotome available in 20mm & 30mm cutting wire options
- Enhanced radio opaque distal tip provides an excellent fluoroscopic visualization
- 5.5Fr distal tip design compatible with 0.035" guidewire




Available In

Sizes: 20mm, 30mm Cutting Wire

Lengths: 180cm, 200cm

Guide Wire Compatibility Upto 0.035"

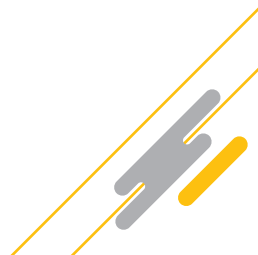

CYTOBRUSH

Cytology Brush

- The brush is designed to get an excellent cytology sample
- The brush is equipped with a metal tip that allows perfect visibility under fluoroscopy



Available In:

Brush Head Diameter: 2mm

Length: 230cm

X-TRACT

Stone Extraction Balloon

- Single balloon provides three distinct diameters which helps to reduce procedure time by eliminating unnecessary exchanges
- Triple lumen catheter provides excellent injection flow while accepting 0.035" guidewire
- Available with injection above and injection below option




Available In:

Injection Port Above, Injection Port Below

Sizes: 9-12-15mm

Length: 200cm


Guidewire Compatibility Upto 0.035"

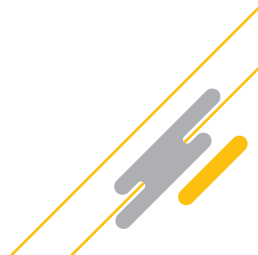

RETRIVO

Stone Extraction Basket

- The gastro stone basket is available in multiple designs and alloy
- Material engineered to maintain the symmetry of the basket after multiple passes
- The unique catheter design provides an excellent manoeuvrability in the biliary anatomy

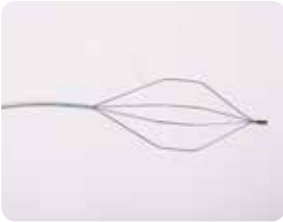

Available In: Hexagonal Shape

Type: Stainless Steel, Nitinol

Open Diameter: 15mm, 20mm, 25mm, 30mm, 35mm

Length: 230cm



Available In: Oval Shape

Type: Stainless Steel, Nitinol

Open Diameter: 15mm, 20mm, 25mm, 30mm, 35mm

Length: 230cm




Available In: Trapezoidal Shape

Type: Stainless Steel, Nitinol

Open Diameter: 25mm, 30mm, 35mm


Length: 230cm

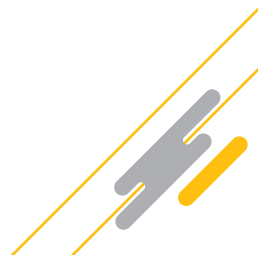

TEFLOFLEX

Plastic Stents

- Used to endoscopically drain bile ducts
- Smooth surface to minimise tissue stress
- Flap & pigtail shape to improve retention
- A rapid flow rate
- Highly radiopaque for better visualization

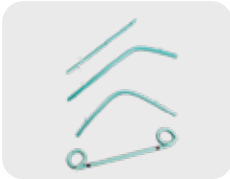

Classification: Standard, Amsterdam, Straight, Single pigtail, Double pigtail, Tannenbaum, Pancreatic with single pigtail, Pancreatic with straight tip.

Material: PTFE

Size: 3-12 Fr

Length: 5-15 cm

LITHOFRAG

Mechanical Lithotripter

- **Variants:** TTS - Through The Scope & OTW - Over the Wire
- **Traction spirals** facilitates controlled, quick & easy lithotripsy
- **Lithotripter baskets** that always breaks at a designed load
- **Reusable & can be steam-sterilized**



Spiral: 3.5 TTS & 4.0mm OTW

Length: 195 cm TTS & 100 cm for TTS

Lithotripter basket size: 25, 30 & 35mm with working length of 400cms

SELF EXPANDING METAL STENTS (SEMS)




EVOLVIA

Metal stents

ESOPHAGEAL SEMS

- **Radial Force**
Braided with single wire to generate uniform radial force and maintain luminal patency
- **Anti-migration mechanism**
Designed with wider proximal and distal flares to reduce chances of migration
- **Prevent tissue in-growth**
Coated with silicone to minimize tissue in-growth
- **Fluoroscopic Visualization**
8 radio opaque markers offering excellent visibility during fluoroscopic visualization
- **Stent placement Accuracy**
Delivery system is reconstrainable up to 95-98% of deployment to aid in repositioning

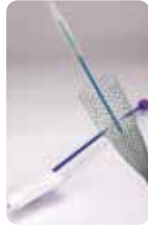

BILIARY SEMS

- **Radial Force**
Braided with single wire for uniform radial force
- **Atraumatic Ends**
Atraumatic ends designed to reduce chances of injury post stent placement
- **Fluoroscopic Visualization**
8 radio opaque markers offering excellent visibility during fluoroscopic visualization
- **Stent Repositioning**
Loop at the proximal end promotes ease in post deployment adjustments
- **Preventing Tissue in growth**
Featured with silicone coating for minimal tissue in-growth
- **Stent Placement Accuracy**
Delivery system is reconstrainable up to 95-98% of deployment to aid in repositioning



ENDOSCOPIC ULTRA SOUND (EUS)



SONOPSY

EUS Guided Needle

FNA Guided needle

- Sharp needle grind helps in precise targeting and sampling
- Excellent needle visibility, as echogenic pattern on the needle extends till the tip of the needle
- Needle constructed with an alloy that helps in superior penetration, excellent pushability and kink resistance
- Needle design engineered for precise cutting capability




FNA (Fine Needle Aspiration)

Size: 19G, 22G, 25G
Length: 141.5cm



FNB (Fine Needle Biopsy)

Size: 19G, 22G, 25G
Length: 141.5cm

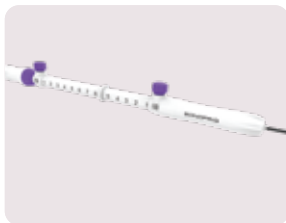

Handle Design



ENTERAL FEEDING



BIOPEG

Percutaneous Endoscopic
Gastrostomy

- **Enhanced grip to prevent migration**
Double O-ring bolster with extra firm grip to prevent migration and facilitate traction removal
- **High strength guidewire**
Smallest & tight knot possible with high strength PU coated guide wire
- **Kink resistant PEG tubing**
Kink resistant medical grade silicone PEG tubing
- Comes with Standard accessories, separate feeding and medication ports pinch clamp to occlude tubing



Available In: Pull Type PEG Kit
Size: 16Fr, 20Fr, 24 Fr, 28Fr
Length: 60cm

STANDARD PEG KIT COMPRISE OF

PULL TYPE PEG TUBE

BOLSTER (WITH O-RING)

ENFIT CONNECTORS

SILICONE FEEDING ADAPTOR

LOOPED GUIDE WIRE

COLD SNARE

SURGICAL FENESTRATED DRAPE

GUIDEWIRE INTRODUCER NEEDLE
WITH SHEATH

PINCH CLAMP

BITE BLOCK


GAUGE SPONGE

SCALPEL BLADE



JEJUNAL FEEDING TUBE

- Used through the channel of the endoscope to provide short duration internal access delivery of nutrition or medication to the small bowel
- Smooth surface for better insertion & less trauma to tissue

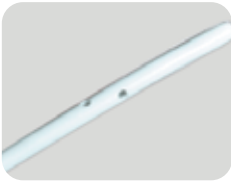

Classification: With Guidewire, Without Guidewire.

Size: 5-16 Fr

Length: 120-200cm

Available as: Normal, PU, Clear



3RD SPACE ENDOSCOPIC ACCESSORIES



TRUSNARE

Polypectomy Snare

- Braided wire design engineered to produce a flexible snare
- Variety of loop configurations available to accommodate a wide range of therapeutic functions
- 360° rotation designed to provide better manoeuvrability in tortuous anatomy



Available In: Hexagonal Shape
Type: Stainless Steel
Open Diameter: 6mm, 13mm, 15mm,
20mm, 25mm, 30mm, 35mm
Length: 230cm



Available In: Oval Shape
Type: Stainless Steel
Open Diameter: 6mm, 13mm, 15mm,
20mm, 25mm, 30mm, 35mm
Length: 230cm



Available In: Crescent Shape
Type: Stainless Steel
Open Diameter: 6mm, 13mm, 15mm,
20mm, 25mm, 30mm, 35mm
Length: 230cm



MISCELLANEOUS



X-PELLER

Loop Basket

- Retrieval net designed to ensure retention of foreign body/sample
- 360° rotation designed facilitate manipulation in tortuous anatomy
- Interwoven net design to maintain the integrity of the net



Available In: Nitinol
Mesh Material: Polystyrene
Open Diameter: 20mm, 30mm, 40mm
Length: 230cm

BIO-NBDC

Nasal Biliary Catheter

- Used for biliary duct drainage through nasal passage and treatment of stricture in the ductile system
- Each set includes a nasal transfer tube and a connecting tube to facilitate irrigation and drainage pigtail shape for better retention



Classification: Straight, Angled with pigtail
Size: 5-12 Fr
Length: 250cm



BIO-BILBAO DOTTER

Enteroclysis Tube

- Used through the channel of the endoscope to provide short duration internal access for delivery of nutrition or medication to the small bowel



Available as: With Guidewire

Size: 8-14 Fr

Length: 150-155cm.



GLOBAL PRESENCE



*Online source



50+ countries
import our quality products

AUSTRALIA

NEW ZEALAND



bioradmedisys.com

Biorad Medisys Pvt. Ltd.

Follow
us on



Corporate Office & Plant:

Survey No.: 48/3 & 48/7, Pashan Sus Rd, Sus Village,

Taluka Mulshi, Pune - 411021. MH, India. Phone: 020-35000200

